

LIMS-Driven Acquisition, Processing, Reporting and Storage of Data from Multiple LC/MS/MS Workstations

David Brant, Jonathan McNally, Joel Usansky, Jim Edwards, Kevin Cook

Thermo Fisher Scientific, San Jose, CA, USA



Overview

Purpose: To develop a LIMS-driven data acquisition and processing software of triple stage quadrupole series instruments.

Methods: A LIMS system that acts as a centralized control point to multiple LC/MS/MS systems in the laboratory has been developed. The LIMS acquires the data and stores it in a database. The data is processed from remote locations and reported. The LIMS controls user access to analytical studies. It also controls which LC/MS/MS instruments can acquire data for the study. The LIMS system is a networked storage location for the data. Remote workstations can integrate and regress the data stored in the LIMS through a web interface, without requiring local installation of LIMS software. This implementation places all analytical study records in a single database simplifying archival.

Results: A LIMS-driven application has been developed to simultaneously handle system security, acquisition, processing, reporting and storage of data from multiple remote triple stage quadrupole series instruments.

Introduction

Historical approaches to interfacing bioanalytical LC/MS/MS instruments with Laboratory Information Management Systems (LIMS) have usually involved the use of two disparate data systems. This has often resulted in the required maintenance of two archival record sets in order to fully document the analytical studies per 21 CFR Part 11 guidelines. This has increased the level of effort required in software validation and direct data analysis, and requires persistent documentation to support other users of the results, such as internal quality assurance offices and external regulatory bodies.

This poster describes a new, integrated approach that results in increased efficiency through the reduction of validation efforts while providing a single archival record set.

Methods

Thermo Watson Bioanalytical LIMS version 7.4 was used as a centralized control point for multiple Thermo Scientific TSQ Quantum and TSQ Vantage systems in the laboratory, the TSQ instruments used the Watson TSQ Module to acquire data directly to the Watson LIMS database. Watson LIMS was used within web browsers from remote locations via the Thermo Web Access Suite to create, review data, process data, and generate reports. Access to analytical studies, and the Watson TSQ Module was controlled by permissions configured within the Watson LIMS. The Watson TSQ Module was also used to send analytical runs back to the Watson Bioanalytical LIMS so the analytical run could be sent to another TSQ series instrument running the Watson TSQ Module.

Results

An application has been developed that controls user access to multiple analytical projects, remote LC/MS/MS instrument workstations, and software functionality. The multiple remote LC/MS/MS systems were added to the LIMS database as separate instrument stacks that could be assigned analytical run sequences. The application allowed analytical run sequences created within the LIMS to be sent to the remote LC/MS/MS systems. The analytical run sequence was assigned and submitted to a particular LC/MS/MS instrument within the LIMS. The analytical run sequences on the remote LC/MS/MS systems were only available to the specific users that were granted access to the analytical run from within the LIMS. These users were then able to submit the analytical run sequence to the particular LC/MS/MS system instrument queue.

FIGURE 1. Configuration of multiple acquisition workstations. Extensible architecture allows multiple workstations to be configured with multiple sample inlets.

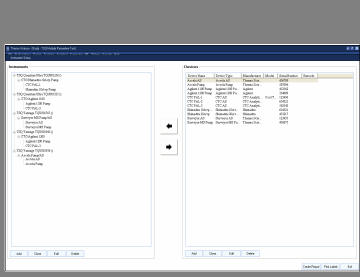


FIGURE 2. Analytical method importation into the database.

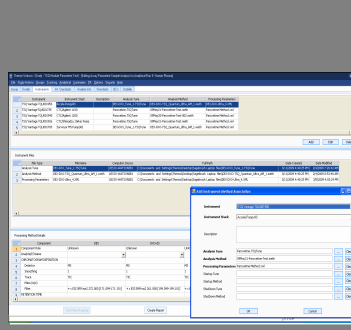


FIGURE 3. Real time display of the currently acquiring sample.



Multiple analytical runs have been successfully submitted to an LC/MS/MS instrument queue and acquired in succession without additional user input. No additional user input is required because each analytical run in the queue begins acquiring data to the LIMS once the preceding analytical run completes its last sample and the raw data is written to the LIMS database.

Analytical run raw data was sent from the remote LC/MS/MS instrument workstations to the LIMS where it was stored, viewed, integrated and then regressed. Viewing and processing of the data was performed from remote workstations through a web interface. This data has then been archived and restored to the database the archive was created from, as well as other separate installations of the LIMS application.

Successful multiple concurrent data acquisitions from many acquisition workstations the LIMS, with users reviewing and processing the acquired data in near real time from multiple concurrent remote locations have been performed. Only those acquisition workstations added to the LIMS database can receive analytical runs and acquire data (Figure 1). Access to both the acquisition workstation module and the remote LIMS application software functionality are managed by the same set of permissions resulting in greater efficiency. Analytical methods used to acquire and process data are developed at an acquisition workstation and then imported into the database (Figure 2). The imported analytical method is associated with a particular instrument at this point, resulting in the correct instrument receiving the analytical run once it is submitted to an acquisition workstation module.

The system provides users with access to both the currently acquiring sample as well as those that have already been acquired. The real time data display of the currently acquiring sample is viewed from the acquisition workstation module at the acquisition workstation (Figure 3), while already acquired samples can be viewed from the remote LIMS application (Figure 4).

FIGURE 4. Review of acquired data and selection of data for reinjection. Selection of samples for reinjection, samples 6 through 9 are selected for reinjection.


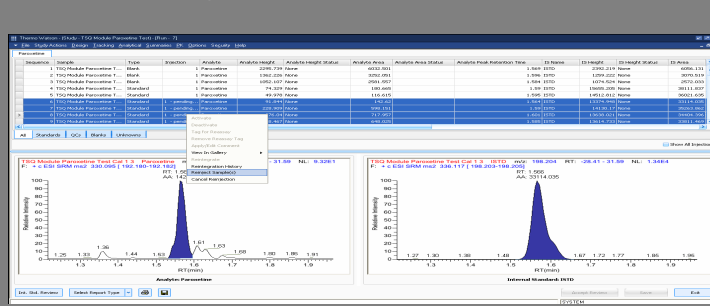



FIGURE 5. Data review of reinjected samples. Original injections are displayed in gray, while the injection used in regression calculations are shown in black.

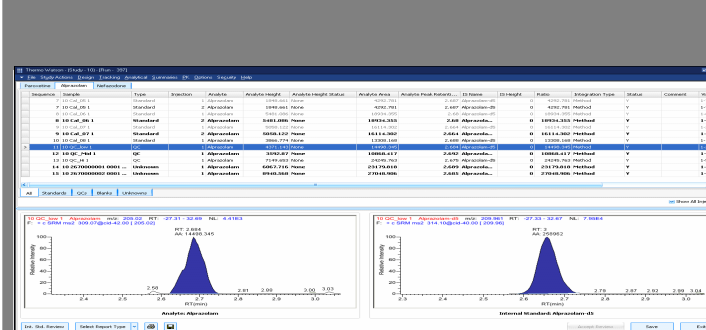

It was important that a sample could be reinjected if required, and then regressed with the other samples in the analytical run without overwriting the original injection. This system is capable of sample reinjection, from both the module at the acquisition workstation if the analytical run is still in progress, or from the remote LIMS processing application once the analytical run is complete. The original injection data is not overwritten and remains in the database where it can be viewed, but cannot be used in the analytical run regression. Multiple reinjections of a sample may be performed, but only the last injection is used in the data regression (Figure 5).

The LIMS application is also capable of reintegrating samples when required. The system allows users to change integration parameters and then apply the new parameters to a single sample or the entire sequence of samples. It is also possible to perform manual integrations, but a manual integration may only be applied to the current sample. The original integration parameters and chromatogram are displayed as a reference in the same window as the modified integration parameters and chromatogram (Figure 6). Integration changes can be displayed and printed in a historical report showing how each integration change affected the chromatogram.


FIGURE 6. Reintegration of acquired data. Original integration parameters and associated chromatogram are displayed on the left of the screen, and the modifications to the chromatogram for the changed integration parameters are displayed on the right side of the screen. Manual integrations are also performed on the right side of the screen, by selecting one of the squares and dragging it to the desired position.

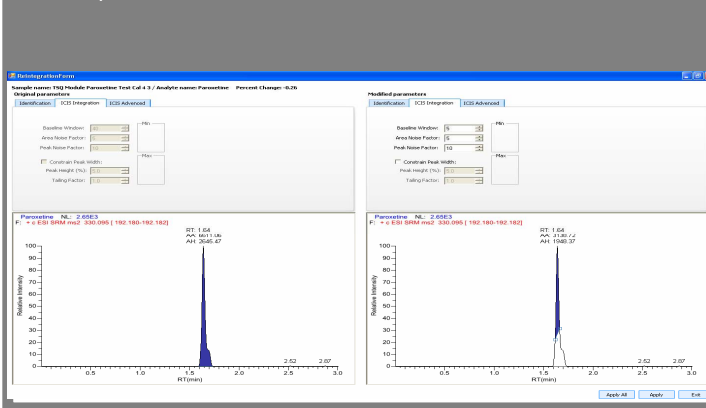

The LIMS application stores the analytical run chromatograms and chromatogram reports within the database. These chromatograms can be viewed using a gallery view function and reports that contain the chromatograms and analytical conditions used to generate the chromatogram can be created (Figure 7). The LIMS application is capable of creating an archive for data retention outside the database. These archives contain all the information regarding the analytical runs, the receiving of samples, raw and processed data, reports, and the individuals that performed each action.

This LIMS-driven application performs functions that have historically been performed by a multiple software packages. Separate software packages have been used for sample handling, data collection, storage, and processing of analytical results. Each application and interface to between applications used in the process required its own validation and maintenance. Previously a change to one application required additional efforts to confirm that the associated existing application interfaces were still functional or resulted in a cascade of changes through multiple applications. This LIMS-driven application approach reduces the number of applications that a laboratory administrator needs to validate and perform administrative functions on resulting in greater efficiency.

FIGURE 7. Review of chromatograms and creation of chromatogram reports.

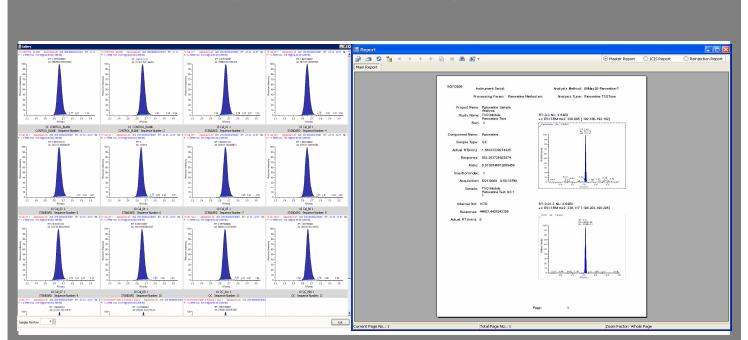

Conclusions

A LIMS-driven application has been successfully developed to simultaneously handle system security, acquisition, processing, reporting and storage of data from multiple remote LC/MS/MS systems. Future work on this project will include the ability to perform analytical method development within the LIMS driven workflow.

All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries.